

OMNIA ACCESSORIES OPERATING THEATRE





MAKE THE RIGHT CHOICE

Accessories such as Instrument Trolley, Stand for Cooling Systems, Surgical Tables and LC Implant Suites, enhance chairside efficiency with standardised procedures, and save time while incrementing the productivity for a complete Chairside Management.

OVER-THE PATIENT SURGICAL TABLE

Omnia over the patient surgical table is made of 18/10 stainless steel, is adjustable in height, and is characterized by a wide and solid supporting surface, which allows the dental surgeon to keep all the material he/she needs to properly work within reach, with no need for additional supports. Furthermore, the dental surgeon may move the working surface over the chest of the patient, so as to avoid the risk of accidental falls of the instruments to the floor. The over the patient surgical table is electropolished by immersion to ensure a lower bacterial adherence. Two out of the 4 wheels are equipped with brakes.



Partcode	Description	Qty
30.E0023.00	Over the patient surgical table with manually adjustable height from 98 to 150 cm (tray: 60 x 90 cm) - 15 kg	1 piece

TABLE ACCESSORY

The removable supporting pole can be mounted onto both Mayo and over-the patient tables and allows positioning the Hand Fuser over the tables.



Partcode	30.E0017.00
Description	Removable supporting pole for Hand Fuser
Qty	1 piece

INSTRUMENT TROLLEY

To complete the equipment of any surgical theatre, Omnia tables are entirely made of stainless steel, with a 4-leg base in die-cast aluminium on castors. Base legs contain heavy alloy hidden weights for a better stability. The tray is large in order to accommodate a great number of surgical instruments. Tables are electropolished by immersion to ensure a lower bacterial adherence.

Partcode	Description	Qty
30.E0026.00	Instrument trolley with edges work top, extra heavy foot anti overturn and 180° rotation, height adjusts from 74 to 117 cm (tray size: 73 x 47 x 2,5 cm) - Weight 18,70 kg. Supporting pole not included.	1 piece



OLEODYNAMIC INSTRUMENT TROLLEY

The oleodynamic Mayo table, which is entirely made of stainless steel, with a 4-leg base in die-cast aluminium on castors, completes the equipment of any dental surgical theatre. Base legs contain heavy alloy hidden weights for a better stability. The oleodynamic column and foot pedal allow a quick adjustment in height of the working surface. The tray, 73 x 47 cm, allows accommodating a large number of surgical instruments. The table is electropolished by immersion to ensure a lower bacterial adherence.

Partcode	Description	Qty
30.E0022.00	Oleodynamic instrument trolley equipped with single tray with raised edges. It is manually adjustable in height from 86 to 120 cm (tray: 73 x 47 x 2,5 cm) - 31,5 kg	1 piece
30.E0018.00	Rotating junction 360° for instrument trolley with oleodynamic pump (optional item)	1 piece
30.E0017.00	Removable supporting pole for Hand Fuser, for positioning the Hand Fuser over Mayo tables and over-the patient tables	1 piece



DETAIL
JUNCTION

STAND FOR COOLING SYSTEMS

The stand with wheeled base is made of stainless steel. It allows positioning and supporting the OMNI-VAC aspiration system, and the wheels allow moving it easily around the surgical room. The tray may be used as a shelf to support the Physiodispenser.



STAND FOR
PHYSIO-
DISPENSER



SUPPORTING
TRAY FOR
HAND FUSER

Partcode	Description	Qty
30.E0014.00	Stand for cooling fluid with tray for Physiodispenser, capacity 10 kg (tray dimensions: 32 x 30 cm) 5 spokes plastic foot	1 piece
30.E0027.00	Removable support for body fluid collection canister that can be added to item 30.E0014	1 piece


30.E0014.00



LC IMPLANT SUITE

The multifunctional surgical trolley is developed to improve the organization of the dental office, it allows to keep the electronic equipment properly stored, to avoid the provisioning of cables and wirings, to improve the level of hygiene in the office, to save time and money in the organization of the work.

LC Implant suite is an essential instrument for those who practices oral surgeries/implants using Phisiodispenser or piezosurgical equipment and for those who practice endodontics and parodontics.

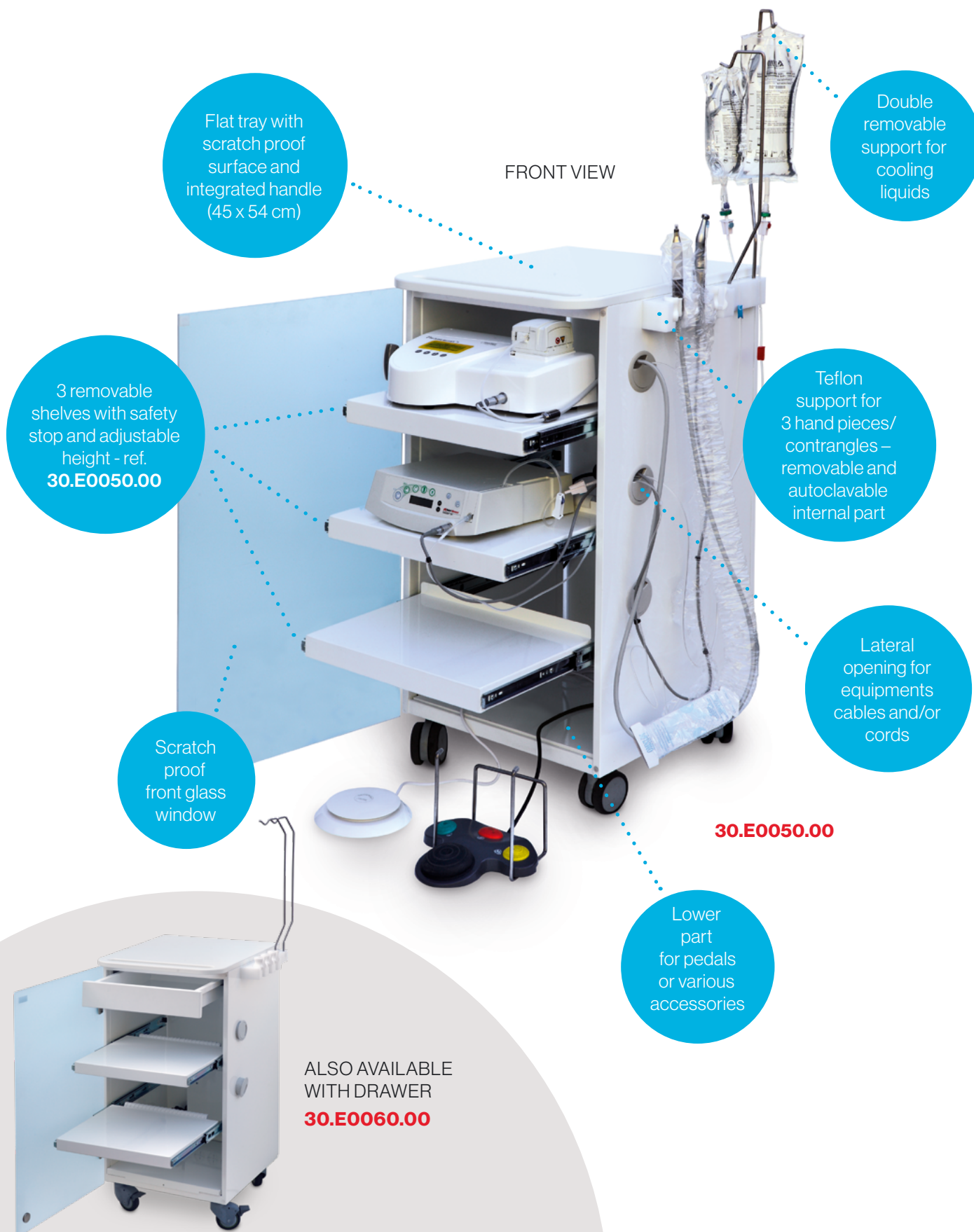



Partcode	Description	Dimensions	Color	Qty
30.E0050.00	LC Implant Suite - 3 removable shelves	50 x 54 x 94 cm	○ white	1 piece
30.E0068.00	LC Implant Suite - 3 removable shelves	50 x 54 x 94 cm	● orange	1 piece
30.E0070.00	LC Implant Suite - 3 removable shelves	50 x 54 x 94 cm	● blue	1 piece
30.E0060.00	LC Implant Suite - 2 removable shelves + 1 drawer	50 x 54 x 94 cm	○ white	1 piece
30.E0072.00	LC Implant Suite - 2 removable shelves + 1 drawer	50 x 54 x 94 cm	● orange	1 piece
30.E0071.00	LC Implant Suite - 2 removable shelves + 1 drawer	50 x 54 x 94 cm	● blue	1 piece



Each LC Implant Suite trolley is equipped with 1 package of tubing extensions for mechanical irrigation systems free of charge (Ref **32.F0100.00** - 10 pcs.)
The use of the extension is necessary for the proper use of the hose set.

For more information on Omnia products or to place an order, call +39 05 24 52 74 53 or visit OmniaSpa.eu

Hu-Friedy Mfg. Co., LLC • European Headquarters • Lyoner Str. 9 • 60528 Frankfurt am Main, Germany • HuFriedyGroup.eu

As some of the products included in this catalogue might not be available in some countries, we invite you to verify the availability of the product with your local authorised dealer or to contact our customer care at the free phone number 00800 48 37 43 39. Available while stocks last. All other company and product names are trademarks of Hu-Friedy Mfg. Co. LLC, its affiliates or related companies, unless otherwise noted.

©2021 Hu-Friedy Mfg. Co., LLC All rights reserved. 90.L9006.00/0921