

# OMNIA DISPOSABLE GOWNS

STERILE AND NON-STERILE PROTECTION



# MAKE THE RIGHT CHOICE



Alongside the spread of more and more diseases and infections in the world, there is an increase in the awareness of all infection control processes that have to be followed during surgical and non-surgical procedures.

We have the solutions that help you minimize the exposure to cross-infection and blood borne pathogens during daily procedures:

## Omnia Gowns

**Disposable gowns are an important step forward in protecting the team. Made of Non-Woven Fabric (NWF), in different grammages and materials depending on their use, they meet all the requirements in terms of protection, comfort and asepsis. Surgical gowns are wrapped in medical paper and contain hand towels.**



## Sterility Guaranteed

The special ethylene oxide (EtO) gas sterilization procedure guarantees a 5-year sterility compared to a few months with the standard autoclave sterilization procedures used in dental practices.

## OUR FABRICS


**Omnia manufactures its line of products using high-performance fabrics.**

Non-Woven Fabric (NWF) materials, which are expressly designed and developed to improve comfort and safety during dental procedures, ensure the utmost hygiene, care and protection to both operators and patients. Their minimum thickness makes the products soft and comfortable and, at the same time, ensures an excellent barrier effect. PE materials are manufactured using special mixtures to ensure safety and hygiene, without impairing the operation of the equipment, where they are used on.

# INNOVATIVE FABRICS

## for improved comfort and safety

## FABRICS USAGE

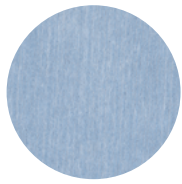

### Sontara®

Thanks to a particular manufacturing process, Sontara® is very soft, drapable and highly breathable. It is the highest qualitative standard ever achieved in medical Non-Woven Fabric technology.




Gowns



Drapes



Caps

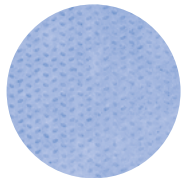

### SMS Microfiber

SMS microfiber (Spunbond + Meltblown + Spunbond) indicates a tri laminate non-woven fabric with two external Spunbond layers and a middle Meltblown layer. SMS microfiber is water-repellent and breathable and guarantees an effective barrier to the passage of fluids and pathogens, yet providing an excellent wearing comfort for the operator.




Gowns



Uniforms



Caps

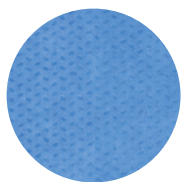


### SMMS Microfiber

The SMMS microfiber is made of spunbonded non-woven fabric, two melt-blown non-woven fabrics and one last spunbonded non-woven fabric that are hot rolled together in the manufacturing process. Water-repellent and breathable as the SMS non-woven fabric, the SMMS Microfiber is comfortable and safe, making these gowns suitable for high-quality isolation performances.



Gowns

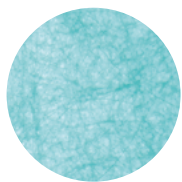

### Spunbond NWF (Non-Woven Fabric)

The low-grammage, soft and light NWF (Non-Woven Fabric) provides real protection for operators and patients. It is used to make gowns for short performances.




Gowns



Caps



Chair Covers

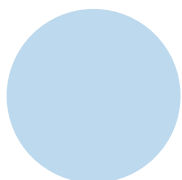

### Non-Toxic Medical Grade PE (polyethylene)


Non-toxic medical grade PE guarantees a complete impermeability to pathogenic and contaminant agents, and makes any procedure safer and easier, without impairing the operation of the equipment used.



Gown Covers



Chair Covers

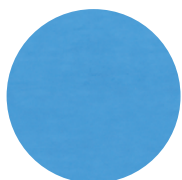


### Absorbent NWF + Liquid proof PE (polyethylene)

The high-grammage Non-Woven Fabric (NWF), melted together with a polyethylene layer, is used to produce surface and patient drapes. It is absorbent and waterproof, tear-proof and comfortable for patients. It ensures a high barrier effect.



Drapes

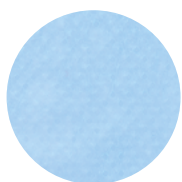


### NWF + Liquid proof PE (polyethylene)

The low-grammage Non-Woven Fabric (NWF), melted together with a liquid proof polyethylene layer, is used to produce surface and patient drapes.



Drapes

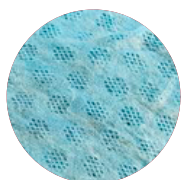

### Aerated NWF (Non-Woven Fabric)

The softness of this NWF (Non-Woven Fabric) and the micro-perforation, which allows ventilation, make it particularly suitable for caps and headwear to be used in dental environments.



Caps


# EXPO MODEL GOWN Sontara®

available **STERILE** and **NON-STERILE**

**Sontara®** is so far one of the most advanced medical non woven fabric available on the market. Its high liquid-repellence, combined with an extraordinary sense of comfort, makes this gown the ideal choice for long working days and the treatment of risk patients. Its main features are: Velcro closure, cotton cuffs, double internal and external side laces, surgically folded and included towels.



**Velcro closure**

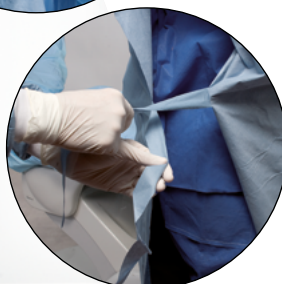

**Double internal and external side laces**



**Cotton cuffs**



**Surgically folded, include 2 hand towels**

22.D1278.00

STERILE	NON-STERILE	DESCRIPTION	PCS
<b>22.D1278.00</b>	<b>30.D1278.00</b>	<b>Expo</b> model gown (70 g) with surgical folding and towel wrapped in medical-grade paper (length: 118 cm) S	12
<b>22.D1275.00</b>	<b>30.D1275.00</b>	<b>Expo</b> model gown (70 g) with surgical folding and towel wrapped in medical-grade paper (length: 130 cm) M/L	12
<b>22.D1276.00</b>	<b>30.D1276.00</b>	<b>Expo</b> model gown (70 g) with surgical folding and towel wrapped in medical-grade paper (length: 133 cm) XL	12
<b>22.D1277.00</b>	<b>30.D1277.00</b>	<b>Expo</b> model gown (70 g) with surgical folding and towel wrapped in medical-grade paper (length: 135 cm) XXL	12



# SPECIAL PLUS MODEL GOWN SMS Microfiber

available **STERILE** and **NON-STERILE**



Water-repellent and breathable, **SMS Microfiber** is used to produce this gown providing an excellent relationship between comfort and performance. Its cotton cuffs, velcro closure, side-fastened waist-belt, surgical folding and internal towels make it the ideal choice for clinical work sessions.

**Velcro closure**



**Cotton cuffs**



**Surgically folded, include 2 hand towels**



22.D2170.00

STERILE	NON-STERILE	DESCRIPTION	PCS
<b>22.D2170.00</b>	<b>30.D2170.00</b>	<b>Special Plus</b> model gown (43 g) with surgical folding and towel wrapped in medical-grade paper (length: 109 cm) S	12
<b>22.D2171.00</b>	<b>30.D2171.00</b>	<b>Special Plus</b> model gown (43 g) with surgical folding and towel wrapped in medical-grade paper (length: 112 cm) M	12
<b>22.D2172.00</b>	<b>30.D2172.00</b>	<b>Special Plus</b> model gown (43 g) with surgical folding and towel wrapped in medical-grade paper (length: 115 cm) L	12
<b>22.D2173.00</b>	<b>30.D2173.00</b>	<b>Special Plus</b> model gown (43 g) with surgical folding and towel wrapped in medical-grade paper (length: 118 cm) XL	12
<b>22.D2174.00</b>	<b>30.D2174.00</b>	<b>Special Plus</b> model gown (43 g) with surgical folding and towel wrapped in medical-grade paper (length: 145 cm) XXL	12

# BLUE SURGICAL GOWN SMMS Microfiber

available **STERILE** and **NON-STERILE**

Water-repellent and breathable as the SMS Non-Woven Fabric, the **SMMS Microfiber** is comfortable and safe, making these gowns suitable for high-quality isolation performances. Its cotton cuffs, velcro closure, side-fastened waist-belt, make it the ideal choice for clinical work sessions.



Velcro closure



Cotton cuffs



Surgically folded

22.D2301.00

STERILE	DESCRIPTION	PCS	NON-STERILE	DESCRIPTION	PCS
22.D2301.00	Surgically folded <b>SMMS</b> gown (35 g) (length: 120 cm) size M	30	30.D2301.00	Essential <b>SMMS</b> gown (35 g) (length: 120 cm) size M	50
22.D2302.00	Surgically folded <b>SMMS</b> gown (35 g) (length: 125 cm) size L	30	30.D2302.00	Essential <b>SMMS</b> gown (35 g) (length: 125 cm) size L	50
22.D2303.00	Surgically folded <b>SMMS</b> gown (35 g) (length: 130 cm) size XL	30	30.D2303.00	Essential <b>SMMS</b> gown (35 g) (length: 130 cm) size XL	50
22.D2304.00	Surgically folded <b>SMMS</b> gown (35 g) (length: 135 cm) size XXL	30	30.D2304.00	Essential <b>SMMS</b> gown (35 g) (length: 135 cm) size XXL	50

Length of gowns may be subject to small variations



# ISOLATION GOWN Spunbond Non-Woven Fabric

available **STERILE** and **NON-STERILE**



**PACK  
OF 50**

Omnia fluid-repellent surgical gowns are the ideal solution for short routine activities, such as the extraction of wisdom teeth, second-stage surgeries and dental scaling, or to supplement the protection and isolation of patient and assistant against pathogens. Available with non-allergic elastic wrist bands or cotton cuffs. (One size only)



**Also available in**

**green**

**22.D1006.00**

**yellow**

**22.D1304.00**



**22.D1302.00**

STERILE	NON-STERILE	DESCRIPTION	DETAILS	COLOR	PCS
<b>22.D1004.00</b>	<b>30.D1004.00</b>	Gown (30 g) with elastic wrist bands (length 120 cm)		green	50
<b>22.D1307.00</b>	<b>30.D1307.00</b>	Gown (30 g) with elastic wrist bands (length 120 cm)		light blue	50
<b>22.D1303.00</b>	<b>30.D1303.00</b>	Gown (30 g) with elastic wrist bands (length 120 cm)		yellow	50
<b>22.D1006.00</b>	<b>30.D1006.00</b>	Gown (30 g) with cotton cuffs (length 120 cm)		green	50
<b>22.D1302.00</b>	<b>30.D1302.00</b>	Gown (30 g) with cotton cuffs (length 120 cm)		light blue	50
<b>22.D1304.00</b>	<b>30.D1304.00</b>	Gown (30 g) with cotton cuffs (length 120 cm)		yellow	50

Length of gowns may be subject to small variations

## GOWN COVER

### Non-Toxic Medical Grade PE & LDPE

available **STERILE** and **NON-STERILE**

The gown cover is perfect for routine practice of all kinds creating high aerosol dispersion levels. The **Non-Toxic Medical Grade PE (polyethylene)** guarantees a complete impermeability to pathogenic and contaminant agents, whereas the **Non-Toxic Medical Grade LDPE (low-density polyethylene)** makes the item waterproof. Both fabrics (PE and LDPE) protect the professional uniform from any accidental splash or contamination.



STERILE	NON-STERILE	DESCRIPTION	COLOR	PCS
<b>22.D1001.00</b>	<b>30.D1001.00</b>	<b>Gown cover</b> (Shell Cover) - LDPE (22,5 g) (length: 105 cm)	light blue	50
<b>22.D1002.00</b>	<b>30.D1002.00</b>	<b>Gown cover</b> (Shell Cover) - PE (22,5 g) (length: 125 cm) size XL, long sleeve with thumb hole	light blue	50 (sterile) 100 (non-sterile)
<b>22.D1009.00</b>	<b>30.D1009.00</b>	<b>Gown cover</b> (Shell Cover) - PE (22,5 g) (length: 125 cm) size XL, long sleeve with thumb hole	light blue	10




## GOWN COVER Spunbond Non-Woven Fabric and UNIFORMS SMS Microfiber

**NON-STERILE** only

**Spunbond NWF gown covers** and **SMS Microfiber uniforms** provide excellent breathability and resistance. Their primarily recommended use is to preserve the gown during daily procedures (gown cover) or as a replacement (professional uniform).



10.D3003.00

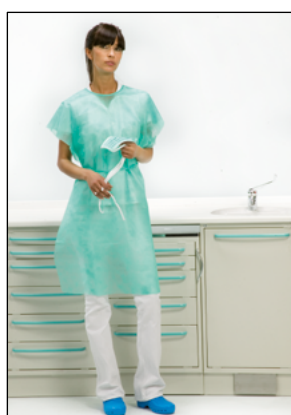

10.D3005.00



30.D2115.00

NON-STERILE	DESCRIPTION	COLOR	PCS
<b>10.D3003.00</b>	<b>Gown cover</b> (length: 110 cm)	white	50
<b>10.D3005.00</b>	<b>Gown cover</b> (length: 110 cm)	green	50
<b>30.D2115.00</b>	<b>Short-sleeve jacket</b> with trousers S	blue	10
<b>30.D2116.00</b>	<b>Short-sleeve jacket</b> with trousers M	blue	10
<b>30.D2117.00</b>	<b>Short-sleeve jacket</b> with trousers L	blue	10
<b>30.D2118.00</b>	<b>Short-sleeve jacket</b> with trousers XL	blue	10



Length of gowns may be subject to small variations

For more information on Omnia products or to place an order, call +39 05 24 52 74 53 or visit [OmniaSpa.eu](https://OmniaSpa.eu)

Hu-Friedy Mfg. Co., LLC • European Headquarters • Lyoner Str. 9 • 60528 Frankfurt am Main, Germany • [HuFriedyGroup.eu](https://HuFriedyGroup.eu)

As some of the products included in this catalogue might not be available in some countries, we invite you to verify the availability of the product with your local authorised dealer or to contact our customer care at the free phone number 00800 48 37 43 39. Available while stocks last. Sontara® is a trademark of Sontara® South Asia. All other company and product names are trademarks of Hu-Friedy Mfg. Co. LLC, its affiliates or related companies, unless otherwise noted.

©2021 Hu-Friedy Mfg. Co., LLC All rights reserved. 90.L9001.00/1121