

# HuFriedyGroup

## Key Opinion Leader

### Dr. Rino Burkhardt

Swiss perio kit developed in collaboration with  
Swiss Perio Group



Focus on:  
**Microsurgery**



Country: **Switzerland**



Favourite Instrument:  
**Microsurgical Anatomical Plier  
SPTPDAPV**



Website:  
**[www.swissperio.com](http://www.swissperio.com)**



Inspired by:  
**Jan Lindhe**



VISIT THE  
**KOL COMMUNITY**



My surgical instruments are the sensory-motor extension of my hands.

Dr. Rino Burkhardt



## SWISS PERIO KIT

### FEATURES AND BENEFITS

#### SPECIFIC DESIGNED INSTRUMENTS


- Less tissue trauma by refined working ends.
- Ultimate manufacturing tolerances for precise handling.
- Easy to clean.

#### BLACK INSTRUMENTS / BLACK LINE EXTENSION

- More contrast between the instrument and the tissue.
- Avoid light reflection.
- More sharpened than standard.

#### BALANCED WEIGHT

More control and precision in specific and delicate microsurgery procedures.

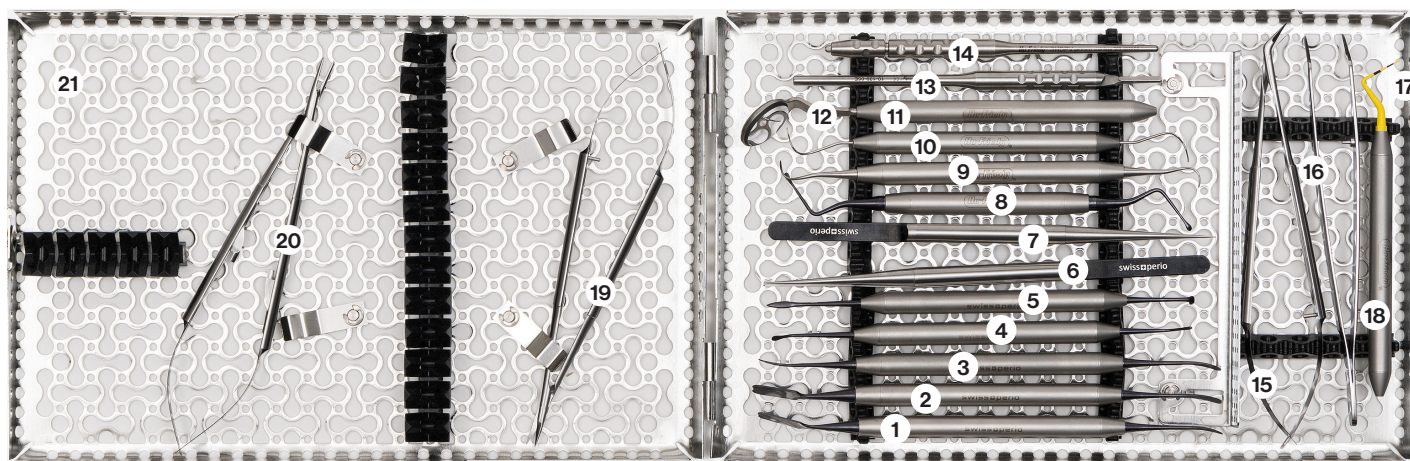

**SPKITDIN**

Complete Swiss Perio Kit

Visit us online at **[HuFriedyGroup.eu](http://HuFriedyGroup.eu)**

SWISS PERIO KIT  
PERIOSTEAL ELEVATORS

BLACK LINE

The **Black Line** is designed to optimise clinical outcomes by delivering efficiency throughout the entire perio and surgical procedure. The Black Line features a performance engineered coating for optimal edge retention, reduced light reflection afforded by the black matte finish and enhanced contrast and visual acuity at the surgical site and underlying tissue.

Three double-ended periosteal elevators are equipped with different working ends based on the variability of oral mucosa, and are designed to cover all kinds of flap elevations from delicate, tiny papilla to coarse palatal masticatory mucosa. The sharpness of the elevators and the ultimate flat, however stiff configuration, guarantee a minimally-traumatic flap mobilization. The blunt working ends with a groove are made to retract full thickness flaps and provide a stable seat on bony surfaces and crests.



**SPP24GSPX**  
Swiss Periosteal #24G  
**Black Line**



**SPP24GSICX**  
Swiss Periosteal #SIC  
**Black Line**



**SPP24GSP2X**  
Swiss Perio Periosteal  
**Black Line**

SWISS PERIO KIT  
PAPILLA ELEVATOR & TISSUE DISSECTORS



**SPPSCHLEEX**  
Swiss Perio Microsurgical  
Plier, straight

The **SPPSCHLEEX** pushes and pulls the papilla. It is indispensable when the interdental area has to be accessed. After intrasulcular sharp dissection of the papilla, its basis and the gingiva of the col area, the interdental mucosa can be elevated and mobilized under ultimate controlled conditions. A key factor for achieving primary wound closure and maintaining wound stability.

In some surgical sites, such as former extraction sockets, the mucosa is firmly adhering to the bone and thus difficult to elevate. The tissue dissector **SPPH2X** allows a controlled dissection of the flap from the underlying hard tissue wound bed. With the opposite end of the instrument, the mucosal tissues can be elevated, directed from posterior to anterior with pulling instead of pushing strokes.

**SPPH2X**  
Swiss Periosteal #PH2  
**Black Line**



SWISS PERIO KIT  
MICROSURGICAL PLIERS



**SPTPDSSPV**  
Swiss Perio Microsurgical  
Plier, straight

The working ends of the surgical forceps **SPTPDSSPV** are configured with micro teeth, allowing to firmly grasp the oral mucosa without traumatizing the flap or graft margins unnecessarily.




**SPTPDAPV**  
Swiss Perio Microsurgical  
Plier, anatomical

The anatomical forceps **SPTPDAPV** are designed for suturing and tightening the knots, but can also be used for manipulating the soft tissue. The very fine and precise configuration of the working ends guarantee a firm seat of the threads and minimize tissue trauma during the surgical procedure.




SWISS PERIO KIT  
MICROSURGICAL SCISSORS

The micro tissue scissors, **SPSPV** with its sharp, slightly curved working end is made for fine resections of soft tissues and trimming of wound margins or mucosal grafts.

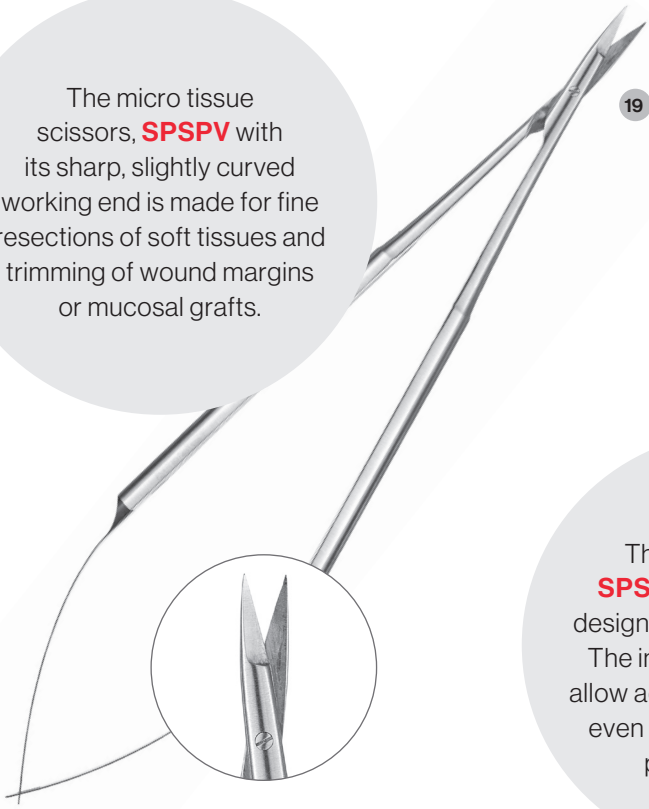

19

The suture scissors **SPSPV45** is exclusively designed for cutting threads. The inclined working ends allow access to difficult areas even far distal to molars in posterior zones.

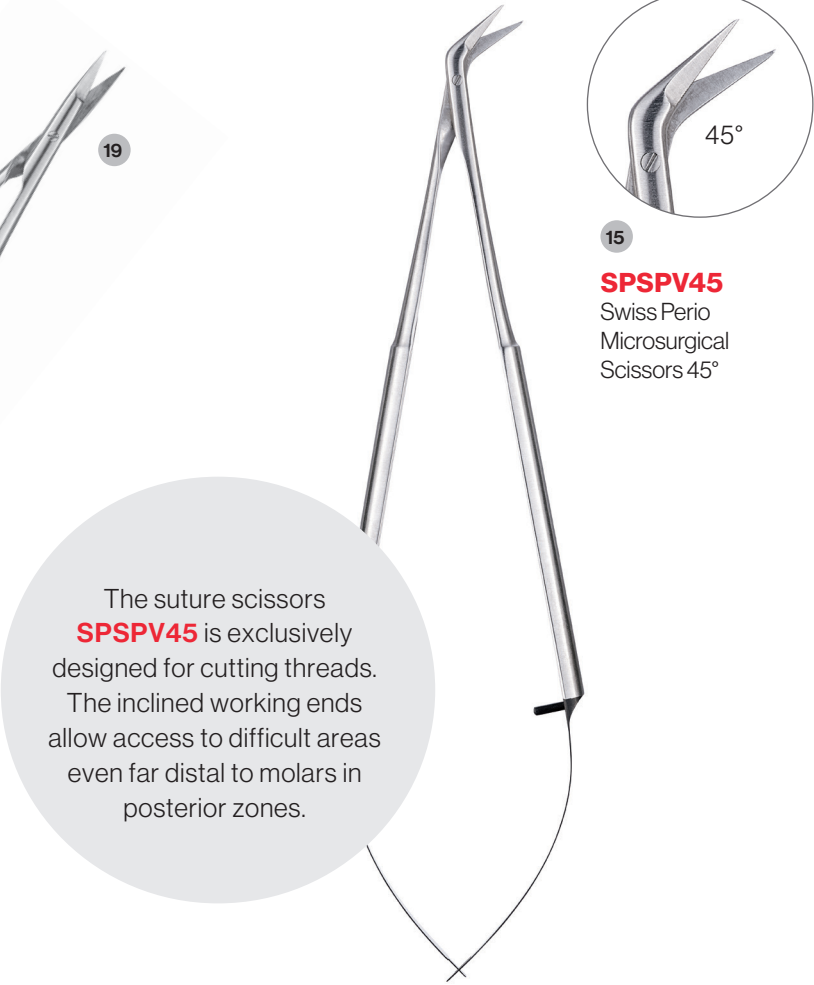


15

**SPSPV45**  
Swiss Perio Microsurgical Scissors 45°

**SPSPV**  
Swiss Perio Microsurgical Scissors

SWISS PERIO KIT  
NEEDLE HOLDER

The needle holder **SPNHDPVN** is completely revised in all details and unmatched in its precision. It allows a firm seat of smaller needles and provides an easy handling of suture diameters in a range of 5-0 to 9-0. The instrument is equipped with a lock that facilitates finger rotation without applying pressure, leading to a more precise execution of finger and hand movements. Additionally, the low locking forces reduce hand tremor while grasping the needles.

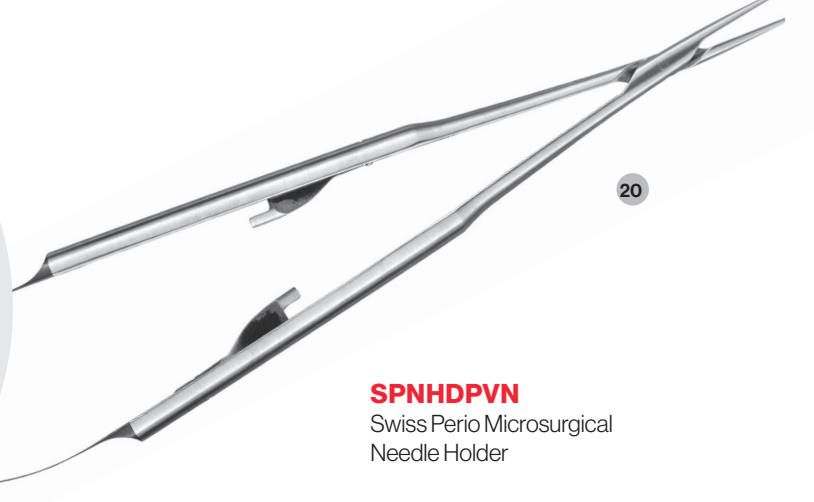

20

**SPNHDPVN**  
Swiss Perio Microsurgical Needle Holder

SWISS PERIO KIT



**DID YOU KNOW?**  
All hinged instruments – forceps, rongeurs, scissors, pliers, hemostats, orthodontic pliers, etc. – should be kept lubricated. Regular use of proper lubricants, like Instrument Lubricant Spray (**ILS**), will prevent rust, corrosion and stiff joints to ensure smooth operation. **All hinged instruments should be sterilised in the open position.**




12

**MIR5DSHDE**  
High Definition Mirror #5 double-sided, European, 1 pc./pack




11

**MHE6CH**  
N.6 Simple Steam Mirror Handle, Smooth Handle




8

**CL84X**  
84 Lucas Double Ended Surgical Curette, **Black Line**



9

**SYG7/8CHE2**  
Curette Younger Good #7/8 hdl #smooth Kotschy



10

**PQ2NCH**  
Q2N Color-Coded Nabers Probe, Smooth Handle

SWISS PERIO KIT



14

**SHDPV**  
Mini Scalpel Blade / Mirror  
Handle Stainless Steel



13

**10-130-05E**  
5 European Style Round  
Scalple Handle




17

**PH6CH**  
Colorvue Probe Handle,  
Smooth Handle Stainless Steel



**PCV12PT**  
12 Colorvue Probe tips,  
12 Replacements Tips

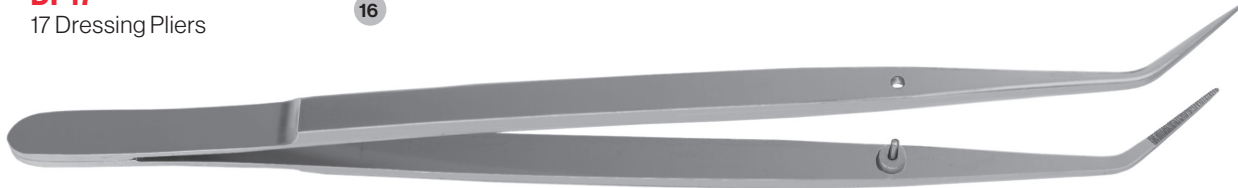

16

**DP17**  
17 Dressing Pliers

SWISS PERIO KIT

CREATE YOUR CUSTOM KIT AND CHOOSE THE CASSETTE WITH YOUR FAVORITE COLOR.

Just add the color coding number to the last part of the partcode (replace “\_” with the number)

|            |   |    |
|------------|---|----|
| Red        | ● | 0  |
| Grey       | ● | 1  |
| White      | ● | 2  |
| Orange     | ● | 3  |
| Purple     | ● | 4  |
| Yellow     | ● | 5  |
| Lavender   | ● | 6  |
| Black      | ● | 7  |
| Blue       | ● | 8  |
| Green      | ● | 9  |
| Ocean Blue | ● | 11 |

WHY YOU SHOULD USE AN INSTRUMENT MANAGEMENT SYSTEM (IMS):

By standardising your instruments with IMS System, additional time and revenue is generated while extending the useful life of your Hu-Friedy instruments. Dental professionals state that by implementing IMS™ Cassettes, they save 5 to 10 minutes per procedure allowing the staff to spend more time on value-added activities\*. In addition, safer transportation of soiled instruments reduces risk of injury and exposure to blood borne pathogens.

\*When compared to single instrument decontamination. Based on market survey results.

DISCOVER  
MORE AT **IMS™**  
WEBPAGE



**IMEDIN18\_S**

IMS™ Infinity Series™ **DIN Surgical Cassette** for 18 instruments, with an accessory area + 4 clips (2 x I1000) for hinged instruments.  
Dimensions: 290 mm x 184 mm x 34 mm



SWISS PERIO KIT WITH CONTAINER CASSETTE:

The Swiss Perio Kit is also available with container cassette:

**SPKITCC**

Complete Swiss Perio Kit with container cassette (2 x **IMECC12\_**)



LEARN MORE  
ABOUT OUR  
**STERILIZATION**  
CONTAINERS.



**IMECC12\_**

IMS™ Infinity Series™ **1/2 DIN Container cassette** with accessory area **12 Instrument**  
Dimensions: 268 x 176 x 34 mm.  
Holds 12 instruments, with handle on lid and filter outlets



## SWISSPERIO EDUCATION



# swiss<sup>+</sup>perio

## EDUCATION

**SwissPerio**, founded by a group of periodontists, is a training institution and think tank offering tailored education in the domain of **periodontics** and **related fields**. All courses and events such as lectures, workshops, discussion groups, clinical clerkships and online activities aim for the prevention and therapy of periodontal and periimplant diseases.

For that reason, **SwissPerio** does not only convey **factual knowledge** and clinical guidelines but critically evaluates research results and translates them into the language of the clinician. As the patient is not a mean nor a median value, a therapy plan cannot be solely based on research findings. Ideally, a treatment should only be applied when both internally and externally valid evidence indicates that it will be useful and beneficial for the particular patient.

Thus, responsible clinical action requires more than just domain knowledge. It requires diagnostic accuracy paired with good medical **decision-making** which is a learned process. Training in critical thinking and decision-making are part of all our courses. These are the tools to link scientific evidence to the clinician's expertise, the patient's individual circumstances and the clinical context and settings.

In clinical practice, success is ultimately dependent on the precise and delicate execution of a given intervention that requires **skillful hands**. Training of manual dexterity and visual-spatial ability are further key aspects in our hands-on courses. Based on representative tasks, clinical procedures are practised stepwise in bench models of different difficulty levels, a learning concept that allows the participants to acquire clinical routine and finally execute the full-blown and technically-demanding interventions safely.

[www.swissperio.com](http://www.swissperio.com)

---

**For more information or to place an order, call at 00800 48 37 43 39 or visit [HuFriedyGroup.eu](http://HuFriedyGroup.eu)**

---

**Free Call:** 00800 48 37 43 39 | **Free Fax:** 00800 48 37 43 40

**E-Mail:** [info@hu-friedy.eu](mailto:info@hu-friedy.eu) | **Website:** [HuFriedyGroup.eu](http://HuFriedyGroup.eu)

   [HuFriedyGroupEurope](https://www.hufriedygroup.eu)

Hu-Friedy Mfg. Co., LLC • European Headquarters • Lyoner Str. 9 • 60528 Frankfurt am Main, Germany • [HuFriedyGroup.eu](http://HuFriedyGroup.eu)  
All company and product names are trademarks of Hu-Friedy Mfg. Co., LLC, its affiliates or related companies, unless otherwise noted.  
©2023 Hu-Friedy Mfg. Co., LLC. All rights reserved. N2310-EN/0323

 **HuFriedyGroup**  
The Best In Practice